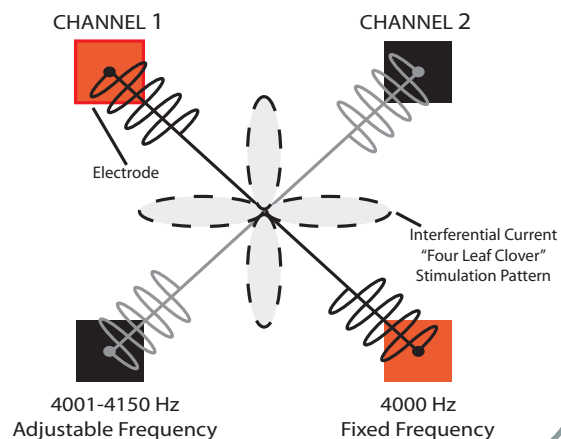


**Elbow Pain**

Quadrupolar placement of the electrodes so the current crosses through the area of pain. Patient should feel the stimulation over the painful region.

**Interferential Quadripolar Diagram****Suggested Treatment Parameters****First Half of Treatment:**

- Beat Frequency: 150 Hz
- Amplitude: Strong but Comfortable (sensory level stimulation)
- Treatment Time: 20 minutes

**Second Half of Treatment:**

- Beat Frequency: 20 Hz
- Amplitude: Maximal Comfort Level (sensory to mild motor level stimulation)
- Treatment Time: 20 minutes

Sessions per day: 1-3

**Suggested Electrode Placement**