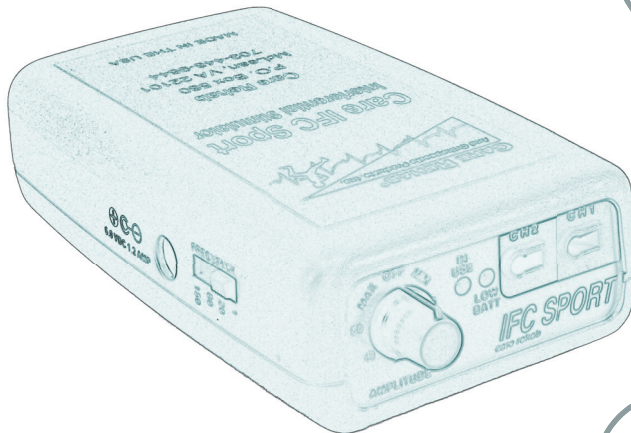
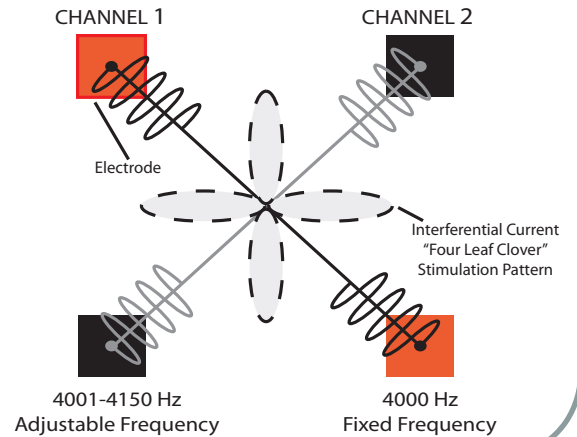
**Neuromas**

Enlarged benign growths of nerves, most commonly between the third and fourth toes. They are caused by tissue rubbing against and irritating the nerves. Pressure from ill-fitting shoes or abnormal bone structure can create the condition as well.

**Interferential Quadripolar Diagram****Suggested Treatment Parameters****First Half of Treatment:**

- Beat Frequency: 80 Hz
- Amplitude: Strong but Comfortable (sensory level stimulation)
- Treatment Time: 20 minutes

**Second Half of Treatment:**

- Beat Frequency: 1 Hz
- Amplitude: Maximal Comfort Level (sensory to mild motor level stimulation)
- Treatment Time: 20 minutes

Sessions per day: 1-3

**Suggested Electrode Placement**