




Bunionectomy

Bunionectomy (or Hallux valgus correction) often results in pain along the medial aspect of the great toe and foot. TENS can be used to decrease pain and increase functional mobility for patients with this condition.

Suggested Electrode Placement

2-pad placement of electrode pads over superior and inferior aspect of great toe

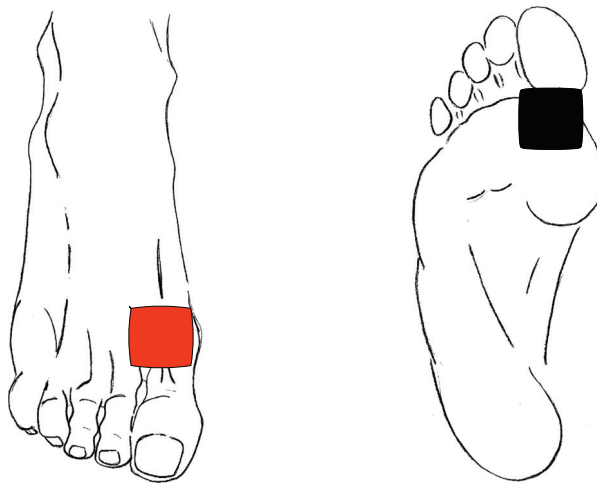

Suggested Treatment Parameters

Mode: Continuous

Pulse rate: 80-100 Hz

Pulse width/duration: 50-100 μ s

Amplitude: Low to moderate level stimulation

Treatment Time: Throughout the day as needed

Mode: Modulated

Pulse rate: 100-120 Hz

Pulse width/duration: 50-150 μ s

Amplitude: Low to moderate level stimulation

Treatment Time: Throughout the day as needed

